

How to explore cancer-related mutations using NGS •



Why is it important to identify mutations in the tumor?

Cancer is primarily caused by genetic alterations, which may be inherited or acquired in somatic cells. These alterations in oncogenes, tumor suppressor genes or DNA repair genes may enable cells to evade growth inhibition in response to growth-inhibitory signals, leading to the development of a tumor. Further mutations may occur in the progeny of the cancer cell, leading to the development of tumor cell invasion and metastasis. The application of NGS provides a high-resolution and global view of the cancer genome in a high-throughput manner, which has revealed a number of novel cancer-related genes and pathways.

Discovery of genetic mutation promotes cancer diagnosis and treatment

A newly identified genetic mutation that occurs frequently in tumor cells may be helpful in cancer diagnosis, monitoring treatment, and prognosis.

Cancer diagnosis and monitoring treatment

Certain mutations in cells can be used for cancer diagnosis, and testing cells for the mutation(s) can be used to check if the cancer is responding to certain treatment. For example, patients with chronic myeloid leukemia (CML) contain a mutated gene BCR-ABL. This mutation is used to confirm the diagnosis.

Drugs targeting genes or gene mutations

Gene mutations may offer drug targets. Drugs can be designed to target the gene, gene mutation, or the protein made by the abnormal gene. For example, patients with HER2-positive breast cancer do not respond to certain chemotherapy drugs. Some drugs such as trastuzumab have been designed to specifically attack HER2-positive cells.

Cancer prognosis

Specific gene mutations can be used to predict which patients tend to have a better or worse outcome, and guide the intensity of treatment. For example, patients whose leukemia cells have a mutation in the FLT3 gene have a worse prognosis than patients without this mutation.

Drugs targeting DNA methylation

Drugs called hypomethylating agents can reverse methylation, which can be used to treat certain cancers in which some genes are abnormally methylated. For example, some genes are found to be methylated in patients with myelodysplastic syndrome. The hypomethylating agents such as decitabine can be useful in treating this disease by decreasing this abnormal methylation.



Discovery of driver genes

One crucial challenge for understanding the molecular mechanisms underlying cancer is to distinguish driver mutations and driver genes to passenger ones. Driver genes contribute to cancer initiation and development, while passenger genes accumulate in cancer cells but do not contribute to the carcinogenesis. Previous researches indicate that driver gene set has two key properties: (i) High coverage in a large number of samples; (ii) High mutual exclusivity (i.e. a single mutation is enough to disturb one pathway). These two properties can be used for the identification of cancer driver genes.

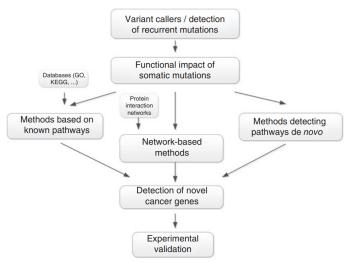


Figure 1. Detection of novel cancer genes.

Table 1. Methods for identification of driver genes in cancer.

Method	Data	Technique
MuSiC	SNV	Statistical test
MutSigCV	SNV	Statistical test
Oncodrive-CLUST	SNV	Clustering
MADGIC	Somatic mutation	Bayesian model
CADD	General variants	Support vector machine
GSEA	Mutation and known pathways	Statistical test
HotNet	Mutation and interaction networks	Heat diffusion process and statistical test
NetBox	Mutation and interaction networks	Community detection
ActiveDriver	Mutations in protein signalling sites	Statistical model
OncodriveFM	Somatic mutations	FM bias
OncodriveFML	Somatic mutations	FM bias
DeepDriver	Somatic mutation	Deep convolutional neural networks



Discovery of pathways

A number of evidences have suggested that pathways usually function cooperatively in cancer initiation and progression. Therefore, exploring the collaboration among different biological pathways and functional modules probably sheds light on the molecular mechanisms underlying carcinogenesis. A number of methods have been developed for the identification of driver gene sets in cancer.


Table 2. Examples of methods for identification of driver pathways in cancer.

Method	Application	
CoMDP	De novo discovers co-occurring mutated driver pathways in cancer by introducing a novel weight function.	
GAMToC	Genomic Alteration Modules using Total Correlation (GAMToC) integrates copy number and mutation data to discover structural module of gene.	
PSMP	Mutual exclusivity-based pairwise search for patterns of somatic mutations.	
Multi-Dendrix	Simultaneously identify multiple driver gene sets in cancer based on the two properties of driver gene sets.	
RME	Detect functional modules based on the patterns of recurrent and mutually exclusive aberrations.	
TiMEx	A generative probabilistic model for de novo detection of mutual exclusive cancer alterations.	
МЕМО	Detects candidate driver subnetworks of aberrant genes with mutually exclusive patterns.	
MEMCover	Uncover pan-cancer dysregulated pathways by combining across tissue type exclusivity with interaction data	
HotNet2	Use insulated network diffusion to detect mutated subnetworks with statistically significant size. Captures the directionality of interactions.	
ComMDP and SpeMDP	De novo identifies cancer common and specific driver gene sets, respectively.	



NGS Approaches

NGS can detect the genomic, transcriptomic, and epigenomic changes including mutations, copy number variations (CNVs), structural variants, fusion genes, differentially expressed genes, and DNA methylation changes, etc.

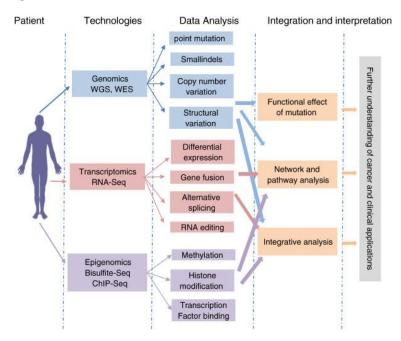


- Cancer Genomics -- Whole genome sequencing, whole exome sequencing, and targeted sequencing (such as gene panel sequencing) can be used to capture genomic changes in tumor tissue. Coupled with powerful bioinformatics, cancer genomics allows researchers to understand oncogenes, tumor suppressor genes, and other genetic factors that contribute to tumor initiation and progression.
- ◆ Cancer transcriptomics -- Monitoring changes in gene expression is important to identify and characterize oncogenes. Related sequencing technologies include RNA-seq, targeted RNA sequencing, single-cell RNA sequencing, and so on. Tools such as WGCNA (weighted correlation network analysis) can be used to construct gene network, allowing to identify driver genes.
- Cancer epigenomics -- Epigenomics technologies (such as whole genome bisulfite sequencing and ChIP-seq) can be used to explore methylation abnormalities and transcription factor binding in tumor, allowing researchers to identify which methylations are crucial for tumor initiation and development.



Our cancer genomics research solutions

CD Genomics offers the sequencing approaches above using next-generation sequencing or longread sequencing technologies to explore multiple types of genetic mutations in tumor, hence understanding the molecular mechanisms underlying carcinogenesis.

- I. The workflow of NGS consists of four steps: library preparation, sequencing, quality analysis, and data analysis and interpretation.
- II. We have advanced sequencers from Roche, Illumina, Ion Torrent, Pacbio, and Oxford Nanopore.
- III. The bioinformatics analyses probably include sequence assembly, function annotation, identification of DEGs (differentially expressed genes), pathway enrichment analysis, pathway association network analysis, discovery of driver genes and pathways, et al.

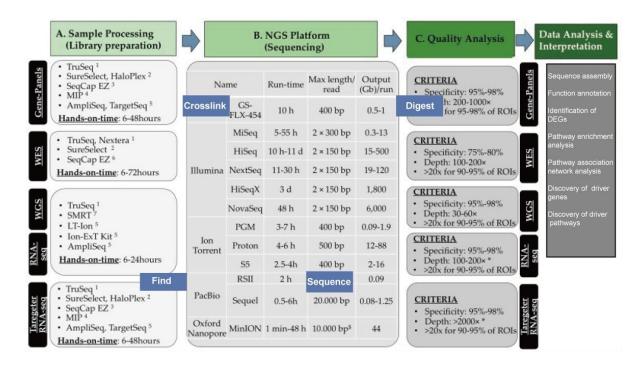

We can tailor this pipeline based on your specific research project.